




Predicta® Bioactive Cement

Predicta Bioactive Cement is designed to fill and seal microgaps with hydroxyapatite formation at the margin. This is said to prevent microleakage, which is otherwise associated with recurrent decay. It is a self-adhesive, self-etching, dual-cure resin cement indicated for luting crowns, bridges, inlays, onlays, and posts (prefabricated metal/nonmetal/fiber posts). It is said to offer outstanding bond strengths across a full range of indications, including zirconia, without the need for printers.

What Makes This Cement Stand Out

By combining all the properties and performance benefits of premium universal resin cements with bioactivity, Predicta Bioactive is said to offer a cement that gives crown margins protection against secondary caries. It also contains 10-Methacryloyloxydecyl dihydrogen phosphate (MDP) for strong bonds to zirconia. Predicta Bioactive Cement is formulated to provide optimal viscosity and consistency for easy handling and easy cleanup. It is highly radiopaque and has a low film thickness for complete seating of restorations with no post-op sensitivity. It is currently available in translucent and A2 shades.

Parkell, Inc
800-243-7446
parkell.com
CIRCLE RS #03

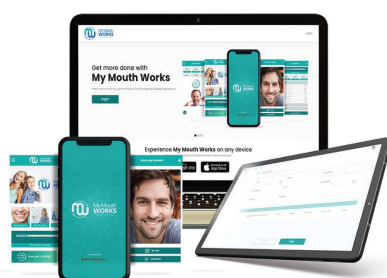

MyMouthWorks

Described as a complete telehealth home care solution, MyMouthWorks is an end-to-end myofunctional therapy program from Innovative Oral Solutions. With the program, clinicians have access to the weekly progress of all their patients including progress reports and images, and clinicians can manage and monitor their patients using a user-friendly web dashboard.

What Sets This Telehealth Platform Apart

With MyMouthWorks, clinicians can conduct an initial patient assessment and then assign 1 of 8 program disorders—open bite, overbite, underbite, crossbite, tongue thrust, temporomandibular joint dysfunction disorders, lingual frenectomy, and sleep apnea. Next, users download the MyMouthWorks mobile app. From there the patient begins a 16-week program and clinicians monitor activity on the Doctors Dashboard and compare weekly progress to baseline. The solution is designed to empower and reward patient adherence, while the doctor-patient notifications and reminders reduce the need for in-office visits. Clinicians can improve patient results with telehealth and experience unlimited ability to treat more patients and increase revenue with low overhead.

Innovative Oral Solutions
760-438-3300
iosmmw.com
CIRCLE RS #04



Bulk EZ Plus™

This new dual-cure, bulk fill composite builds on Bulk EZ's enhanced esthetics and polishability with new features. It is an easy-to-place dual cure composite that combines unlimited cure depth with flowable cavity adaptation, high strength, and enhanced durability. Bulk EZ PLUS is now said to be suitable for a wider range of clinical situations through its thixotropic properties. It reaches a higher degree of conversion, offers better polish retention, and is stackable. Bulk EZ PLUS is available in a variety of shades.

What Makes This Composite Different

Utilizing the same Intellitek Technology for unlimited cure depth and microleakage prevention as the original, Bulk EZ PLUS now incorporates 80-nm spherical zirconia fillers. This is said to enhance the esthetics with better shade-matching and polishability, ensuring greater wear resistance.

Zest Dental
800-262-2310
zestdental.com
CIRCLE RS #05