

MucoHardTM+

Hard Denture Reline Material




Hard Reliner with Low Exothermic Temperature for Improved Patient Comfort.

MucoHard+ is a self-cured, pink gingival-colored, hard relining material used for in-office, durable long-term relines, repairs, and for additions to the borders of full or partial acrylic dentures. Automix cartridge dispensing promotes an ideal mix every time for convenient chairside hard relines. The viscosity of MucoHard+ enables easy placement and spreading, complete seating of the denture base back onto the upper or lower arch, and detailed capture of the topography to ensure a snug fit.

MucoHard+ is a biocompatible, tasteless, odorless, MMA-free hard reline material for full and partial dentures. It has a low exothermic temperature for improved patient comfort during the initial intraoral self-curing stage, and when set, a color consistent with the shade of the most common lab processed denture base acrylics.

Saves time for you and your patient, while saving the denture from a trip to the lab.

- ✓ Low exothermic temperature for improved patient comfort
- ✓ Efficient and easy-to-use
- ✓ Eliminates powder/liquid mixing
- ✓ Auto-polymerizing resin formula sets without a curing light
- ✓ Pink polymer blends beautifully with most denture base shades

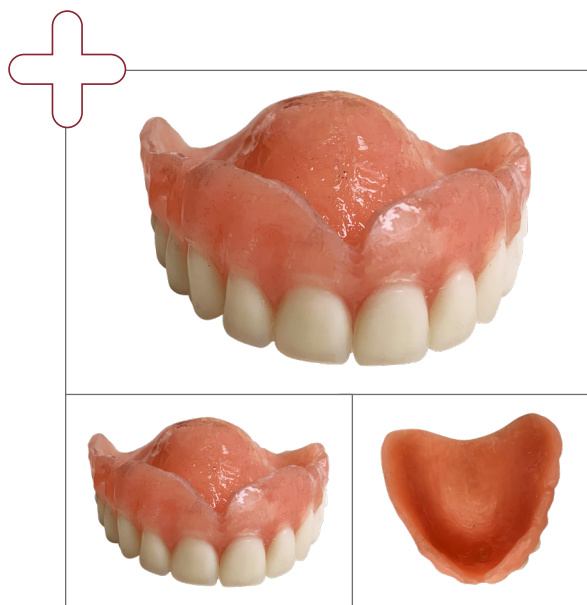

A MucoHard+ in-office hard denture reline: Used to restore fit, function and comfort, in less than 45 minutes.

SKU# DESCRIPTION

S640	MucoHard+ Hard Denture Reline Kit; Includes (1) 50ml cartridge of auto-polymerizing resin, (1) 15ml bonding liner, (6) Mixing tips, brushes and pipettes
------	---

ACCESSORIES & REFILLS

S641	MucoHard+ Reline Resin Cartridge (50ml)
------	---

S642	MucoHard+ Bonding Liner (15ml)
------	--------------------------------

S496	Blue-Green Base Mixing Tip (30 tips)
------	--------------------------------------



For more information, visit Parkell.com or contact your authorized Parkell distributor.

parkell[®]